


Series 2000XL

2-sectioned versatile treatment tables with electric motorized height adjustment ... for examination, treatment, therapy ...

Model 2006-00XLE and 2007-00XLE white powder coated (RAL 9010)
 Model 2006-03XLE and 2007-03XLE white aluminium powder coated (RAL 9006)
 Model 2006-04XLE and 2007-04XLE silver aluminium powder coated (RAL 9007)

► Model 2806-..XLE corresponds with model 2006-..XLE, but with nose opening

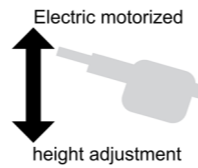

All pictures:
 Model 2006-00XLE
 with option: Movable on
 comfort twin castors,
 antistatic (1186-10);
 upholstery colour:
 64-174 sky

Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Head part can be hinged vertical
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive (art. no. 1999, mandatory equipment)
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless up- (+30°) and downwards (-90°) adjustable by gas spring
- Length head part: 450 mm
- Fixed foot part
- Length foot part: 1500 mm
- Upholstery thickness: 60 mm



Technical data:

- Electronical height adjustment: 520 mm to 960 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: Model 2006-..XLE: 700 mm | Model 2007-..XLE: 800 mm
- Length of lying surface: min. 1550 mm - max. 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
 Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
 Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
 Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
 Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Additional motor control options



Art. no. 1920E
 Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
 Foot switch (low voltage)



Art. no. 1940E
 1 foot switch fixation mounted at underframe



Art. no. 1950E
 Foot switch rail, lateral, on both sides, security adjustment

Paper roll holder



Art. no. 1270K / 1280K
 Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



Art. no. 1270F / 1280F
 Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

Side guards



Art. no. 1080
 1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 1060 mm



Art. no. 1040
 1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm

Wall protection bumpers



Art. no. 1175 / 1176
 4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Nose opening



Art. no. 1610
 Nose opening, in head part

Art. no. 1620
 Upholstery part to close

Upholstery options

Art. no. 1770 / 1780
 Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
 Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm