

Series 2000XL


3-sectioned versatile treatment tables with electric motorized height adjustment
... for examination, treatment, therapy ...

Model 2008-00XLE and 2009-00XLE white powder coated (RAL 9010)
Model 2008-03XLE and 2009-03XLE white aluminium powder coated (RAL 9006)
Model 2008-04XLE and 2009-04XLE silver aluminium powder coated (RAL 9007)

► **Model 2808-..XLE corresponds with model 2008-..XLE, but with nose opening**



Both pictures:
Model 2008-00XLE
with options:
Movable on twin castors with
central lock (1180N), nose
opening and upholstery part to
close (1610, 1620);
upholstery colour:
64-701 stone blue

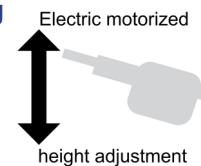

Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Foot part can be hinged vertical
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive (art. no. 1999, mandatory equipment)
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+50°) and downwards (-45°) adjustable by gas spring
- Length head part: 450 mm
- Fixed middle part
- Length middle part: 1290 mm
- Foot part stepless up- (+30°) and downwards (-90°) adjustable by gas spring
- Length foot part: 450 mm
- Upholstery thickness: 60 mm



Technical data:

- Electronical height adjustment: 520 mm to 960 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: Model 2008-..XLE: 700 mm | Model 2009-..XLE: 800 mm
- Length of lying surface: min. 1790 mm - max. 2190 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology




Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin
castors Ø 100 mm,
each with total lock,
basic height -15 mm

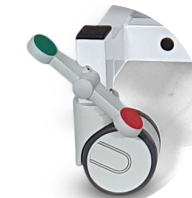

Art. no. 1180N
Movable on twin castors,
Ø 100 mm with central
lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin
castors, Ø 125 mm with
central lock,
basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Additional motor control options



Art. no. 1920E
Hand switch (low voltage),
with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at
underframe



Art. no. 1950E
Foot switch rail, lateral, on both sides,
security adjustment

Paper roll holder



Art. no. 1270K/F / 1280K/F
Paper roll holder, mounted at foot
part (not protruding),
for width 700 mm / 800 mm

Wall protection bumpers



Art. no. 1175 / 1176
4 wall protection bumpers Ø 150 mm, inclu-
ding mounting,
for width 700 mm / 800 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close

Upholstery options

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for
width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for
width 700 mm / 800 mm, thickness 70 mm