


Series 2000XL

3-sectioned versatile treatment tables with with electric motorized height adjustment
... for examination, treatment, therapy ...

Model 2010-00XLE and 2011-00XLE white powder coated (RAL 9010)
Model 2010-03XLE and 2011-03XLE white aluminium powder coated (RAL 9006)
Model 2010-04XLE and 2011-04XLE silver aluminium powder coated (RAL 9007)

► Model 2807-..XLE corresponds with model 2010-..XLE, but with nose opening

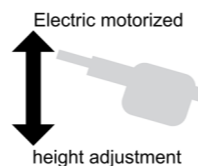

Model 2010-00XLE
with options:
Movable on comfort twin castors
with central lock (1185);
upholstery colour:
64-717 ddeep sea

Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Foot part can be hinged vertical
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive (art. no. 1999, mandatory equipment)
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless adjustable from 0° to +60° by gas spring
- Length head part: 450 mm
- Fixed middle part
- Length middle part: 1.050 mm
- Foot part stepless up- (+30°) and downwards (-90°) adjustable by gas spring
- Length foot part: 450 mm
- Upholstery thickness: 60 mm



Model 2010-00XLE with options:
Movable on comfort twin castors with central lock (1185), paper roll holder (1270K/F)

Technical data:


- Electronical height adjustment: 530 mm to 970 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: Model 2010-..XLE: 700 mm | Model 2011-..XLE: 800 mm
- GesamtLength of lying surface: min. 1550 mm - max. 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai@ Pandoria Plus with velvet-smooth surface
- Skai@ Toronto EN with innovative staynu-technology

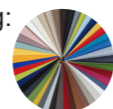

Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Additional motor control options



Art. no. 1920E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at underframe



Art. no. 1950E
Foot switch rail, lateral, on both sides, security adjustment

Paper roll holder



Art. no. 1270F/K / 1280F/K
Paper roll holder, mounted under head part (70 mm protruding), for width 700 mm / 800 mm



Art. no. 1270K/F / 1280K/F
Paper roll holder, mounted at foot part (not protruding), for width 700 mm / 800 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close

Wall protection bumpers



Art. no. 1175 / 1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Side guards



Art. no. 1080-60
1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 820 mm

Upholstery options

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm