

# Series 2000XL

2-sectioned versatile treatment tables with electric motorized height adjustment  
... for examination, treatment, therapy ...

Model 2052-00XLE and 2053-00XLE white powder coated (RAL 9010)  
Model 2052-03XLE and 2053-03XLE white aluminium powder coated (RAL 9006)  
Model 2052-04XLE and 2053-04XLE silver aluminium powder coated (RAL 9007)

► Model 2910-..XLE corresponds with model 2052-..XLE, but with nose opening



3 upper pictures:  
Model 2052-00XLE  
with options:  
Movable on comfort twin  
castors with central lock  
(1185); upholstery colour:  
64-716 black

## Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Lying surface inclinable to 22° (Trendelenburg) by electromotive drive in 18 sec.
- Operation of electric motor to adjust lying surface by 1 hand switch
- Operation of electric motor for height adjustment by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drives (art. no. 1999, mandatory equipment)
- Electric motor for height adjustment: Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



## Upholstery:

Lying surface 2-sectioned:

- Head part stepless up- (+45°) and downwards (-20°) adjustable by gas spring
- Length head part: 500 mm, optional with nose opening and upholstery part
- Fixed foot part
- Length foot part: 1.500 mm
- Upholstery thickness: 60 mm



Model 2052-00XLE with options: Movable on twin castors with central lock (1180N), 2 side guards lateral lowerable (1080); upholstery colour: 64-174 sky

## Technical data:


- Electronical height adjustment: 540 mm to 980 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Lying surface inclinable to 22° (Trendelenburg) by electromotive drive in 18 sec.
- Width of lying surface: Model 2052-..XLE: 700 mm | Model 2053-..XLE: 800 mm
- Length of lying surface: 2000 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



## Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

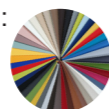

Available as an option (see pricelist for additional options)

## Mobility



Art. no. 1120N-10  
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N  
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185  
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N  
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



## Additional motor control options



Art. no. 1920E  
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E  
Foot switch (low voltage)



Art. no. 1940E  
1 foot switch fixation mounted at underframe



Art. no. 1950E  
Foot switch rail, lateral, on both sides, security adjustment

## Paper roll holder



Art. no. 1270K / 1280K  
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



Art. no. 1270F / 1280F  
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

## Nose opening



Art. no. 1610  
Nose opening, in head part

Art. no. 1620  
Upholstery part to close

## Side guards



Art. no. 1080  
1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 1060 mm

## Upholstery options

Art. no. 1770 / 1780  
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782  
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm