


Series 2000XL

2-sectioned versatile treatment tables with hydraulic height adjustment
... for examination, treatment, therapy ...

Model 2052-00XL/H and 2053-00XL/H white powder coated (RAL 9010)
Model 2052-03XL/H and 2053-03XL/H white aluminium powder coated (RAL 9006)
Model 2052-04XL/H and 2053-04XL/H silver aluminium powder coated (RAL 9007)

► Model 2810-..XL/H corresponds with model 2052-..XL/H, but with nose opening

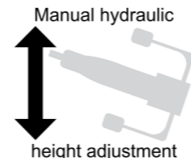

All pictures:
Model 2052-04XL/H
with option: Movable on comfort twin castors with central lock (1185) paper roll holder (1270K), 2 side guards lateral lowerable (1080) (picture shows chromium-plated push handle, basic push handle is RAL 9006 powder coated); upholstery colour: 64-717 deep sea

Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Lying surface inclinable to 12° (Trendelenburg) by gas spring
- Push handle at foot part facilitates the operation of manual inclination, RAL 9006 powder coated
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment

Upholstery:

Lying surface 2-sectioned:

- Head part stepless up- (+45°) and downwards (-35°) adjustable by gas spring
- Length head part: 500 mm
- Fixed foot part
- Length foot part: 1500 mm
- Upholstery thickness: 60 mm





Technical data:

- Manual foot operated height adjustment: 550 mm to 970 mm by a hydraulic pump with a pressure capacity of 10000 N
- Lying surface inclinable to 12° (Trendelenburg) by gas spring
- Width of lying surface: Model 2052-..XL/H: 700 mm | Model 2053-..XL/H: 800 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load

Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

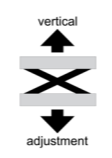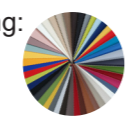

Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Paper roll holder



Art. no. 1270K / 1280K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm

Wall protection bumpers



Art. no. 1175 / 1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Nose opening



Art. no. 1610
Nose opening, in head part
Art. no. 1620
Upholstery part to close

Side guards



Art. no. 1080
1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 1060 mm

Upholstery options

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for width 700 mm/ 800 mm, thickness 80 mm
Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm