

Series 2000XL



4-sectioned multipurpose table with electric motorized height adjustment

... for examination, rectoscopy, gynaecology ...

Model 2640-00XLE white powder coated (RAL 9010)

Model 2640-03XLE white aluminium powder coated (RAL 9006)

Model 2640-04XLE grey aluminium powder coated (RAL 9007)

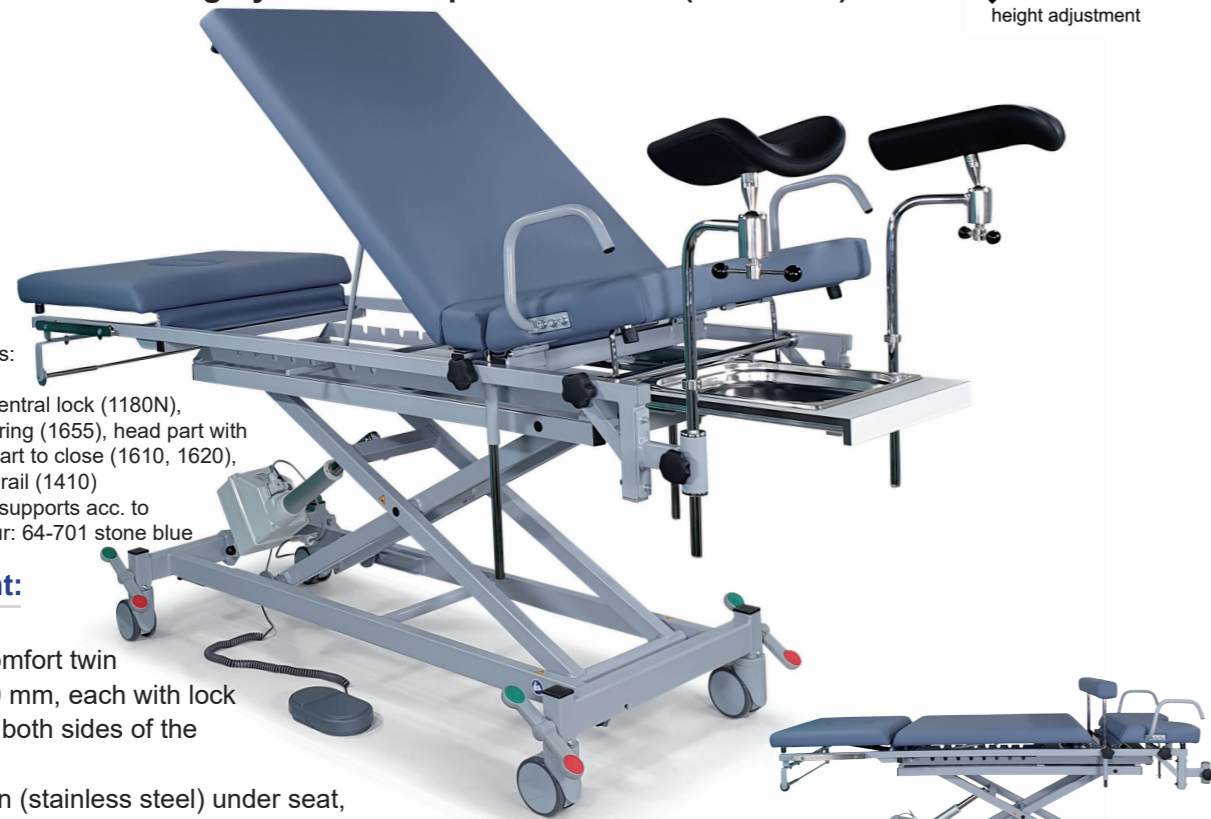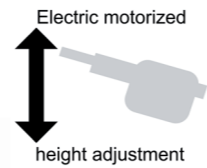

All pictures:
Model 2640-03XLE with options:
loss of basic mobility (1182),
movable on twin castors with central lock (1180N),
head part adjustable by gas spring (1655), head part with
nose opening and upholstery part to close (1610, 1620),
paper roll holder (1265K), side rail (1410)
upper picture additionally: Leg supports acc. to
Goepel (1520); Upholstery colour: 64-701 stone blue

Basic equipment:

- Movable on comfort twin castors, Ø 100 mm, each with lock
- Hand grips on both sides of the seat section
- Drawable basin (stainless steel) under seat, basin support: HPL panel
- Including clamps for leg supports (leg supports = available as an option)
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive (art. no. 1999, mandatory equipment)
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment

Upholstery:

Lying surface 4-sectioned:

- Head part 30° positive adjustable by 2 metal rachets
- Length head part: 550 mm
- Middle part 50° positive adjustable by rack system
- Length middle part: 930 mm
- Rectoscopy insert to lift up and fixable, width: 100 mm
- Seat section with half-round cut-out, manually inclinable (20°)
- Length seat section: 400 mm
- Thickness of upholstery: 80 mm



Technical data:

- Electronical height adjustment: 620 mm to 1060 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 650 mm
- Length of lying surface: 1980 mm
- Load capacity: 220 kg; according to IEC 60601 tested with fourfold static load

kg

Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 1182
Without basic mobility



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close

Additional motor control options



Art. no. 1920E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at underframe



Art. no. 1950E
Foot switch rail, lateral, on both sides, security adjustment



Paper roll holder

Art. no. 1265K
Paper roll holder, mounted under head part (not protruding), for with of lying surface 700 mm

Side guards



Art. no. 1040
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm



Art. no. 1005
1 side guard to fit on side rail, with attachment clamp, chromium-plated, length: 800 mm



Head part options

Art. no. 1655
Head part stepless adjustable by gas spring from 0° to +60°

Add-on elements



Art. no. 1410
1 piece of side rail, length: 350 mm



Art. no. 1430
1 attachment clamp for side rail, suitable for pipe Ø 18 mm



Art. no. 1450
1 Infusion pole for side rail, without attachment clamp

Leg supports



Art. no. 1510
1 pair of leg supports



Art. no. 1520
1 pair of leg supports according to Goepel

Upholstery options

Art. no. 1770
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm

Art. no. 1772
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm

