

# Series 2000XL


## 4-sectioned multipurpose table with hydraulic height adjustment

... for examination, rectoscopy, gynaecology ...

Model 2640-00XL/H white powder coated (RAL 9010)

Model 2640-03XL/H white aluminium powder coated (RAL 9006)

Model 2640-04XL/H grey aluminium powder coated (RAL 9007)

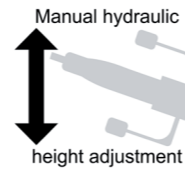

2 upper pictures:  
Model 2640-00XL/H  
with option:  
Paper roll holder (1265K),  
side rail (1410),  
leg supports (1510);  
upholstery colour: 64-174 sky

(Model -01, top frame  
chromium-plated  
no longer available)



### Basic equipment:

- Movable on comfort twin castors, Ø 100 mm, each with lock
- Hand grips on both sides of the seat section
- Drawable basin (stainless steel) under seat, basin support: HPL panel
- Including clamps for leg supports (leg supports = available as an option)
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 4-sectioned:

- Head part 30° positive adjustable by 2 metal rachets
- Length head part: 550 mm
- Middle part 50° positive adjustable by rack system
- Length middle part: 930 mm
- Rectoscopy insert to lift up and fixable, width: 100 mm
- Seat section with half-round cut-out, manually inclinable (20°)
- Length seat section: 400 mm
- Thickness of upholstery: 80 mm



### Technical data:

- Manual foot operated height adjustment: 640 mm to 1060 mm by a hydraulic pump with a pressure capacity of 10000 N
- Width of lying surface: 650 mm
- Length of lying surface: 1980 mm
- Load capacity: 220 kg; according to IEC 60601 tested with fourfold static load

### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

### Mobility



Art. no. 1182  
Without basic mobility



Art. no. 1180N  
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm

### Nose opening



Art. no. 1610  
Nose opening, in head part

Art. no. 1620  
Upholstery part to close

### Paper roll holder



Art. no. 1265K  
Paper roll holder, mounted under head part (not protruding), for with of lying surface 700 mm

### Side guards



Art. no. 1040  
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm



Art. no. 1005  
1 side guard to fit on side rail, with attachment clamp, chromium-plated, length: 800 mm

### Head part options



Art. no. 1655  
Head part stepless adjustable by gas spring from 0° to +60°

### Add-on elements



Art. no. 1410  
1 piece of side rail, length: 350 mm



Art. no. 1430  
1 attachment clamp for side rail, suitable for pipe Ø 18 mm



Art. no. 1450  
1 Infusion pole for side rail, without attachment clamp

### Leg supports



Art. no. 1510  
1 pair of leg supports



Art. no. 1520  
1 pair of leg supports according to Goepel

### Upholstery options

Art. no. 1770  
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm

Art. no. 1772  
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm

