

Tables for echocardiography


3-sectioned table for echocardiography with electric motorized height adjustment

... for examination, treatment ...

Model 2690-00XLE white powder coated (RAL 9010)

Model 2690-03XLE white aluminium powder coated (RAL 9006)

Model 2690-04XLE grey aluminium powder coated (RAL 9007)

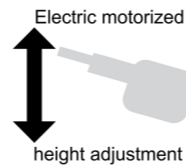

All pictures:
Model 2690-04XLE
with option:
Movable on twin castors with
central lock 1187; with
special feature: foot bar;
upholstery colour
64-175 silver grey



Basic equipment:

- Table movable on twin castors, Ø 75 mm, each with lock
- In middle part: Cut-out for echocardiography and upholstery part to close (in patient supine position left)
- Middle part stepless positive and negative adjustable by electromotive drive, control by hand switch
- Possible positions including: Roof position, Trendelenburg, Reverse Trendelenburg
- Plug connector for potential equalization in acc. to DIN 42801
- Electromotive drive control of height adjustment by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment

Upholstery:

Lying surface 3-sectioned:

- Head part stepless positive (+45°) and negative (-30°) adjustable by gas spring
- Length head part: 400 mm, with nose opening and upholstery part to close
- Middle part (length: 580 mm) stepless positive (+15°) and negative (-15°) adjustable by electromotive drive
- Middle part with cut-out for echocardiography and upholstery part to close, measures: L. 300 mm x W. 250 mm
- Foot part stepless positive (+62°) and negative (-12°) adjustable by gas spring
- Length foot part: 970 mm
- Upholstery thickness: 60 mm



Technical data:

- Electronical height adjustment: 630 mm to 1070 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 220 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 1182
Without basic mobility



Art. no. 1187
Movable on twin castors with central lock, Ø 100 mm, lower underframe, table height 540 mm – 980 mm

Additional motor control options



Art. no. 1920E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at underframe

Nose opening



Art. no. 1615-EF
Without nose opening

Wall protection bumpers



Art. no. 1175
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm

Paper roll holder



Art. no. 1270F/K
Paper roll holder, mounted under head part (not protruding), for width 700 mm



Art. no. 1270F
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm

Add-on elements



Art. no. 1410
1 piece of side rail, length: 350 mm



Art. no. 1005
1 side guard to fit on side rail, with attachment clamp, chromium-plated, length: 800 mm

Upholstery options

Art. no. 1772
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm

