


Heavy-weight tables

2-sectioned heavy-weight table with electric motorized height adjustment

... for examination, treatment, therapy ...

- Model 7209-00E white powder coated (RAL 9010)
- Model 7209-03E white aluminium powder coated (RAL 9006)
- Model 7209-04E grey aluminium powder coated (RAL 9007)

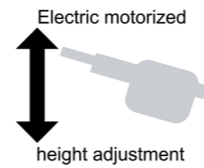

Model 7209-00E with options:
Movable on heavy-weight twin castors with central lock (1189), paper roll holder (1280K); Upholstery colour: 64-701 stone blue

Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Extremely high load capacity: 400 kg
- Electric motorized height adjustment by special heavy-weight AC-motor
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive (art. no. 1999, mandatory equipment)
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless adjustable by gas spring from 0° to +40°
- Length head part: 500 mm
- Fixed foot part
- Length foot part: 1450 mm
- Thickness of upholstery: 80 mm



Model 7209-00E with options:
Movable on heavy-weight twin castors with central lock (1189), foldable side guards (1008)




Technical data:

- Electric motorized height adjustment: 500 mm to 850 mm in 36 sec. by special heavy-weight AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 800 mm
- Length of lying surface: 1950 mm
- Load capacity: 400 kg; according to IEC 60601 tested with fourfold static load

Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

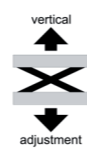

Available as an option (see pricelist for additional options)

Mobility



Art. no. 1189
Movable on heavy-weight twin castors with central lock, Ø 125 mm

Additional motor control options



Art. no. 1920E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at underframe

Paper roll holder



Art. no. 1280K
Paper roll holder, mounted under head part (not protruding), for width of lying surface 800 mm



Art. no. 1280F
Paper roll holder, mounted at foot part (not protruding), for width of lying surface 800 mm

Side guards



Art. no. 1080
1 side guard foldable, chromium-plated, completely made of steel tube, length: 1060 mm

Wall protection bumpers



Art. no. 1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 800 mm


Head part options



Art. no. 1610
Nose opening

Art. no. 1620
Upholstery part to close

Storage basket

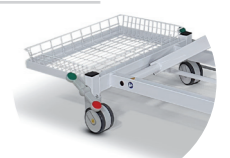

Art. no. 1850
Storage basket with mounting, measures:
L. 530 x W. 370 x H. 80 mm

Upholstery options

Art. no. 1782
Upholstery medium-soft (cold foam cover), for width 800 mm, thickness 70 mm