


Heavy-weight tables

3-sectioned heavy-weight table with electric motorized height adjustment

... for examination, treatment, therapy ...

- Model 7255-00E white powder coated (RAL 9010)
- Model 7255-03E white aluminium powder coated (RAL 9006)
- Model 7255-04E grey aluminium powder coated (RAL 9007)

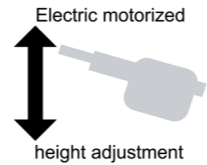

All pictures:
Model 7255-00E with options:
Movable on heavy-weight twin
castors with central lock (1189),
Paper roll holder (1280K);
Upholstery colour: 64-708 cherry

Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Extremely high load capacity: 400 kg
- Electric motorized height adjustment by special heavy-weight AC-motor
- Electromotive drive control by foot switch (option: hand switch)
- Head part and foot part positive adjustable by electromotive drive, controlled by hand switch
- Positioning of an adipose patient is easily possible due to electric motorized adjustment of lying surface segments and height
- Integrated safety function for automatic blocking of electromotive drives (art. no. 1999, mandatory equipment)
- Height adjustment by electromotive drive with freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment

Upholstery:

Lying surface 3-sectioned:

- Head part stepless adjustable from 0° to +38° by electro motive drive
- Length head part: 550 mm
- Fixed middle part
- Length middle part: 410 mm
- Foot part stepless adjustable from 0° to +30° by electro motive drive
- Length foot part: 1030 mm
- Thickness of upholstery: 80 mm




Technical data:

- Electric motorized height adjustment: 540 mm to 890 mm in 36 sec. by special heavy-weight AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 800 mm
- Length of lying surface: 1990 mm
- Load capacity: 350 kg; according to IEC 60601 tested with fourfold static load

Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

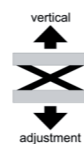

Available as an option (see pricelist for additional options)

Mobility



Art. no. 1189
Movable on heavy-weight twin castors with central lock, Ø 125 mm

Additional motor control options



Art. no. 1920E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at underframe

Wall protection bumpers



Art. no. 1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 800 mm

Paper roll holder



Art. no. 1280K
Paper roll holder, mounted under head part (not protruding), for with of lying surface 800 mm



Art. no. 1280F
Paper roll holder, mounted at foot part (not protruding), for width of lying surface 800 mm

Head part options



Art. no. 1610
Nose opening

Art. no. 1620
Upholstery part to cose

Storage basket



Art. no. 1850
Storage basket with mounting, measures:
L. 530 x W. 370 x H. 80 mm

Upholstery options

Art. no. 1782
Upholstery medium-soft (cold foam cover), for width 800 mm, thickness 70 mm