


Series 4800E

3-sectioned therapy table with electric motorized height adjustment

- Model 4815-00E white powder coated (RAL 9010)
- Model 4815-03E white aluminium powder coated (RAL 9006)
- Model 4815-04E grey aluminium powder coated (RAL 9007)

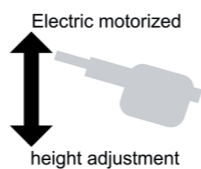

All pictures:
Model 4815-00E
with options:
Movable on twin castors with central lock (046), surrounding foot control (095E), paper roll holder (010F); upholstery colour: 64-708 cherry

Basic equipment:

- Lying surface 3-sectioned
- Foot part stepless adjustable (+75°) to sitting position
- Table not movable, including height adjustment screws to compensate for uneven ground
- High adjustment range
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive (art. no. 1999, mandatory equipment)
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment

Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+50°) and downwards (-40°) adjustable by gas spring
- Length head part: 450 mm; with nose opening and upholstery part to close
- Middle part: fixed
- Length middle part: 530 mm
- Foot part stepless upwards adjustable (+75°) to sitting position by gas spring
- Length foot part: 1020 mm
- Thickness of upholstery: 60 mm

Technical data:


- Electronical height adjustment: 450 mm to 1000 mm in 18 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 2000 mm
- Load capacity: 200 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

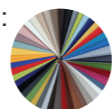

Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 046
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering

Additional motor control options



Art. no. 090E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 091E
Foot switch (low voltage)



Art. no. 092E
1 foot switch fixation mounted at underframe



Art. no. 094E
Foot switch rail, lateral, on both sides, security adjustment:
Press the switch rail down = the table is lifted
Lift the switch rail up = the table is lowered



Art. no. 095E
Foot switch rail, surround, security adjustment:
Press the switch rail down = the table is lifted;
Lift the switch rail up = the table is lowered

Paper roll holder



Art. no. 010K
Paper roll holder, mounted under head part (not protruding), for width 700 mm, not in case of 3-sectioned head part (Art. no. 056)



Art. no. 010F
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm



Art. no. 055
Head part 3-sectioned, arm rests downwards adjustable from -25° to +45° and removable (art. no. 010K not possible) measures:
Head part: L. 500 mm x W. 370 mm
Arm rests: L. 420 mm x W. 160 mm

Arm rests



Art. no. 088
1 arm rest, fixed, swivelling and removable, additional revolving and tiltable

Nose opening



Art. no. 1615-EF
Loss of nose opening and upholstery part to close

Upholstery options

Art. no. 007
Width of upholstery 800 mm

Art. no. B1770
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm

Art. no. B1772
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm

