





TABLES FOR ECHOCARDIOGRAPHY FOR EXAMINATION AND TREATMENT

2023.1

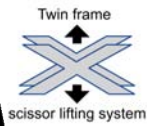

TABLES FOR SURGERY AND HOSPITAL

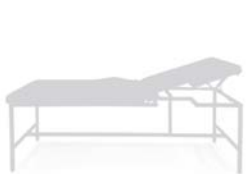

...WITH FIXED HEIGHT



... WITH HEIGHT ADJUSTMENT

Tables for echocardiography



Durable ✓



Robust ✓



Easy-care ✓




Cut-out for echocardiography and self fixing upholstery part ✓




High load capacity ✓

Tables with height adjustment additionally:

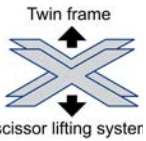

Twin frame scissor lifting system ✓

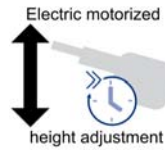

High adjustment speed ✓



Proven Dewert-quality ✓



Content

Good reasons for a table for echocardiography 3

Tables for echocardiography with fixed height

2-sectioned tables for echocardiography:

- Model 117 / 177 6

Tables for echocardiography with height adjustment

Upholstery part(s) positive adjustable

2-sectioned tables for echocardiography:

- Model 2550XLE / 2551XLE 8
- Model 2550XL/H / 2551XL/H 10

3-sectioned tables for echocardiography:

- NEW** • Model 2552XLE 12
- Model 2552XL/H 14

Tables for echocardiography with height adjustment

Upholstery part(s) positive and negative adjustable

- Model 2690XLE 16
- Model 2695XLE 18

Tables for echocardiography with height adjustment and tilt function

- Model 2095XLE 20



All Models with **cut-out for echocardiography**:

- Cut-out for echocardiography and self fixing upholstery part to close
- Conical shape prevents the unintended slipping out of upholstery part




Tables for echocardiography

Good reasons for a table for echocardiography




Models with fixed height

- Frame made of square tube, 30 x 30 x 2 mm
- With cross brace between the feet, longitudinal brace and central support for better stability
- Load capacity: 250 kg

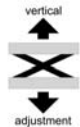


Electric motorized height adjustment

- Electric motorized height adjustment by AC-motor with a pressure capacity of 10000 N and constant adjustment speed even in case of high load
- EPR-safety power cable with strain relief, length: 4 m
- **Integrated safety devices (according to the newest advice of BfArM)**
 - Automatic blocking of electromotive drive
 - Active freewheeling downwards, automatic declutch in case of resistance
- Plug connector for potential equalization in acc. to DIN 42801
- Load capacity: 250 kg

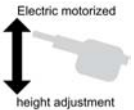

Hydraulic height adjustment

- Manual foot operated height adjustment by a hydraulic pump with a pressure capacity of 10000 N
- Pedals for foot-controlled height positioning on both sides of the table
- Load capacity: 250 kg

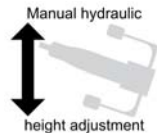

Dewert-hygiene-equipment

Enables optimal cleaning and disinfection

- Artificial leather underneath the upholstery
- Vent holes underneath the upholstery, closed by special caps with valve function
- Easy hygiene because of an open construction: good accessibility of components, no plastic coverings which are breakable in the ageing process
- Option: Innovative artificial leather Skai® Toronto EN with staynu



Tables for echocardiography

Good reasons for a table for echocardiography ...



Upholstery:

- Upholstery thickness: 60 mm, Model 117 / 177: 80 mm
- Multiple glued, high-strength laminated wood panel (thickness: 18 mm) no breaking chipboards in case of pressure on the edges
- Two-layer homogenous polyurethane foam, density 46 and 75, Ökotex-certified
- 2 different variants for artificial leather eligible for upholstery covering:
 - Skai® Pandoria Plus with velvet-smooth surface
 - Skai® Toronto EN with innovative staynu-technology



Variants:

... with fixed height and 2-sectioned lying surface:



- Head part stepless upwards (+55°) adjustable by gas spring
- Longitudinal brace for additional stability
- Galvanised steel hinges connect the upholstery parts
- Frame welded as one unit (powder coated) or composed of two units (chromium-plated)

... with height adjustment:



Technology:



- Extremely stable **twin frame scissor lifting system**, with high torsional stiffness due to thick-walled square tubes with big dimensions, also for under and upper frame and scissor

- High stability due to massive 20 mm axes




- Widely-spaced feet and a high deadweight guarantee stability

- Robust frame construction which is proven in demanding everyday life of hospitals for years



- Abrasion-resistant guide system at under and upper frame due to the use of roller bearings

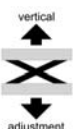



Tables for echocardiography

Good reasons for a table for echocardiography ...

... with height adjustment and 2-sectioned lying surface:



- Head part with cut-out for echocardiography positive stepless adjustable to +50° by gas spring with tubular steel frame for high torsion rigidity

... with height adjustment and 3-sectioned lying surface:

- Head part with cut-out for echocardiography positive stepless adjustable to +50° by gas spring with tubular steel frame for high torsion rigidity
- Fixed middle part
- Foot part positive adjustable by gas spring
- Pressure stable gas spring ensures working at foot part **without vibrations**
- Operating lever to manipulate the gas spring easily accessible



... with height adjustment and 3 adjustable upholstery sections:

- Head part steples positive and negative adjustable by gas spring
- Middle part stepless positive and negative adjustable by electromotive drive
→ Electromotive drive control by hand switch
- With one or two cut-out(s) for echocardiography and upholstery part(s) to close, depending on model



- Foot part stepless positive and negative adjustable:
Possible positions including:
Roof position, Trendelenburg,
Reverse Trendelenburg



... with height adjustment and tilt function:

- Head part stepless positive adjustable by gas spring
- Independent electric motorized height adjustment and tilt function
→ Electromotive drive control of tilt function by hand switch

Tables for echocardiography

2-sectioned table for echocardiography with fixed height ... for examination, treatment ...

Model 117-00 and 177-00 white powder coated (RAL 9010)

Model 117-01 and 177-01 chromium-plated

Model 117-03 and 177-03 white aluminium powder coated (RAL 9006)

Model 117-04 and 177-04 grey aluminium powder coated (RAL 9007)



All pictures:
Model 177-04
in basic equipment;
upholstery colour 64-153
platin



Basic equipment:
Head part adjustment by
gas spring



Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Cut-out for echocardiography and upholstery part to close (in patient supine position left)
- Possibility to pin on a paper roll holder on both traverse sides
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless upwards (+55°) adjustable by gas spring
- Length head part: 800 mm
- Cut-out for echocardiography and upholstery part to close, in patient supine position left, measures: L. 300 mm x W. 250 mm
- Fixed foot part
- Length foot part: 1150 mm
- Upholstery thickness: 80 mm



Technical data:


- Length of lying surface: 1950 mm
- Width of lying surface: Model 117-...: 650 mm | Model 177-...: 800 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

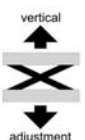

Available as an option (see pricelist for additional options)



Mobility



Art. no. 1120-10
Movable on 4 castors,
Ø 100 mm, each with total
lock



Art. no. 1140-10
Movable on 4 castors,
Ø 125 mm, each with
total lock



Art. no. 1180N
Movable on twin castors,
Ø 100 mm with central lock
(not in case of chromium-
plated frame)



Art. no. 1185
Movable on twin castors,
Ø 125 mm with central lock
(not in case of chromium-plated
frame)

Paper roll holder



Art. no. 1210 / 1220
Paper roll holder, for width
650 mm / 800 mm

Upholstery options

Art. no. 1765 / 1780
Soft covering for pleasant resting
comfort, for width 650 mm / 800 mm,
thickness 80 mm

Art. no. 1762 / 1782
Upholstery medium-soft (cold foam
cover), for width 650 mm / 800 mm,
thickness 70 mm

Side guards



Art. no. 1080-50
1 side guard lateral lowerable,
chromium-plated, completely made of steel
tube, length: 1060 mm,
length on cut-out side: 820 mm



Art. no. 1040
1 side guard tiltable, chromium-plated,
completely made of steel tube,
length: 700 mm



Art. no. 1005
1 side guard to fit on side rail, with
attachment clamp, chromium-plated,
length: 800 mm

Add-on elements



Art. no. 1410
1 piece of side rail,
length: 350 mm

Art. no. 1420
1 piece of side rail,
length: 1300 mm



Art. no. 1430
1 attachment clamp for
side rail, suitable for pipe:
Ø 18 mm



Art. no. 1440
1 arm rest for side rail,
without attachment clamp



Art. no. 1450
Infusion pole, without
attachment clamp



Art. no. 1850
Storage basket with mounting,
measures:
L. 530 x W. 370 x H. 80 mm


Tables for echocardiography

2-sectioned table for echocardiography with electric motorized height adjustment ... for examination, treatment ...

Model 2550-00XLE and 2551-00XLE white powder coated (RAL 9010)

Model 2550-03XLE and 2551-03XLE white aluminium powder coated (RAL 9006)

Model 2550-04XLE and 2551-04XLE grey aluminium powder coated (RAL 9007)

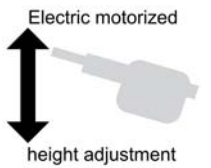

2 upper pictures:
Model 2550-00XLE
with options: Movable on twin
castors with central
lock (1180N),
paper roll holder (1270K);
upholstery colour
64-717 deep sea



Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- In head part: Cut-out for echocardiography and upholstery part to close (in patient supine position left)
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Plug connector for potential equalization in acc. to DIN 42801
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless upwards (+50°) adjustable by gas spring
- Length head part: 800 mm
- In head part: Cut-out for echocardiography and upholstery part to close, in patient supine position left, measures: L. 300 mm x W. 250 mm
- Fixed foot part
- Length foot part: 1150 mm
- Upholstery thickness: 60 mm



Model 2550-00XLE with options:
Movable on comfort twin castors, each with lock (1120N-10), storage basket (1850), electric motorized head part adjustment (1665-10)

Technical data:


- Electronical height adjustment: 500 mm to 940 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: Model 2550-...XLE: 700 mm | Model 2551-...XLE: 800 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandora Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

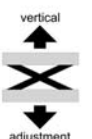

Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering

Additional motor control options



Art. no. 1920E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at underframe

Wall protection bumpers



Art. no. 1950E
Foot switch rail, lateral, on both sides, security adjustment



Art. no. 1170
1 wall protection bumper Ø 150 mm

Head part options



Art. no. 1665-10
Electric motorized head part adjustment to +60°

Paper roll holder



Art. no. 1270K / 1280K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



Art. no. 1270F / 1280F
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

Upholstery options

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

Add-on elements



Art. no. 1410
1 piece of side rail, length: 350 mm



Art. no. 1850
Storage basket

Side guards



Art. no. 1080-50
1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 1060 mm, length on cut-out side: 820 mm



Art. no. 1040
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm



Art. no. 1005
1 side guard to fit on side rail, with attachment clamp, chromium-plated, length: 800 mm


Tables for echocardiography

2-sectioned table for echocardiography with hydraulic height adjustment ... for examination, treatment ...

Model 2550-00XL/H and 2551-00XL/H white powder coated (RAL 9010)

Model 2550-03XL/H and 2551-03XL/H white aluminium powder coated (RAL 9006)

Model 2550-04XL/H and 2551-04XL/H grey aluminium powder coated (RAL 9007)

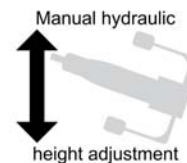

All pictures:
Model 2550-00XL/H
with options: Movable on twin castors with central lock (1180N), paper roll holder (1270K); upholstery colour 64-708 cherry



Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- In head part: Cut-out for echocardiography and upholstery part to close (in patient supine position left)
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless upwards (+50°) adjustable by gas spring
- Length head part: 800 mm
- In head part: Cut-out for echocardiography and upholstery part to close, in patient supine position left, measures: L. 300 mm x W. 250 mm
- Fixed foot part
- Length foot part: 1150 mm
- Upholstery thickness: 60 mm

Technical data:

- Manual foot operated height adjustment: 520 mm to 940 mm by hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: Model 2550-..XL/H: 700 mm | Model 2551-..XL/H: 800 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Paper roll holder



Art. no. 1270K / 1280K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



Art. no. 1270F / 1280F
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

Add-on elements



Art. no. 1410
1 piece of side rail, length: 350 mm



Art. no. 1850
Storage basket

Side guards



Art. no. 1080-50
1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 1060 mm, length on cut-out side: 820 mm



Art. no. 1040
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm



Art. no. 1005
1 side guard to fit on side rail, with attachment clamp, chromium-plated, length: 800 mm

Wall protection bumpers



Art. no. 1170
1 wall protection bumper Ø 150 mm

Upholstery options

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm



Tables for echocardiography

3-sectioned table for echocardiography with electric motorized height adjustment ... for examination, treatment ...

Model 2552-00XLE white powder coated (RAL 9010)

Model 2552-03XLE white aluminium powder coated (RAL 9006)

Model 2552-04XLE grey aluminium powder coated (RAL 9007)

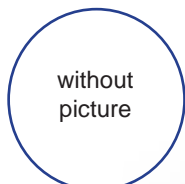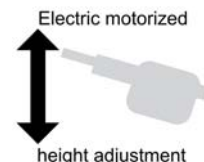

Model 2552-00XL/H with options:
Movable on twin castors with central lock (1180N),
paper roll holder (1270K);
upholstery colour 64-711 ash grey



Frame in case of electric motorized height adjustment



Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- In head part: Cut-out for echocardiography and upholstery part to close (in patient supine position left)
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Plug connector for potential equalization in acc. to DIN 42801
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless upwards (+50°) adjustable by gas spring
- Length head part: 800 mm
- In head part: Cut-out for echocardiography and upholstery part to close, in patient supine position left, measures: L. 300 mm x W. 250 mm
- Fixed middle part, length: 300 mm
- Foot part stepless upwards (+30°) adjustable by gas spring
- Length foot part: 850 mm
- Upholstery thickness: 60 mm



Technical data:

- Electronical height adjustment: 540 mm to 980 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm, **optional** 800 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandora Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering

Additional motor control options



Art. no. 1920E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at underframe

Wall protection bumpers



Art. no. 1950E
Foot switch rail, lateral, on both sides, security adjustment



Art. no. 1170
1 wall protection bumper Ø 150 mm

Head part options



Art. no. 1665-10
Electric motorized head part adjustment to +60°

Paper roll holder



Art. no. 1270K / 1280K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



Art. no. 1270F / 1280F
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

Upholstery options

Art. no. 2552-800
Width of upholstery: 800 mm

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

Add-on elements



Art. no. 1410
1 piece of side rail, length: 350 mm



Art. no. 1850
Storage basket

Side guards



Art. no. 1080-50
1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 1060 mm, length on cut-out side: 820 mm



Art. no. 1040
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm



Art. no. 1005
1 side guard to fit on side rail, with attachment clamp, chromium-plated, length: 800 mm


Tables for echocardiography

3-sectioned table for echocardiography with hydraulic height adjustment ... for examination, treatment ...

Model 2552-00XL/H white powder coated (RAL 9010)

Model 2552-03XL/H white aluminium powder coated (RAL 9006)

Model 2552-04XL/H grey aluminium powder coated (RAL 9007)

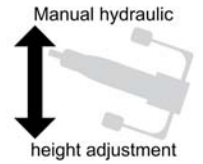

All pictures:
Model 2552-00XL/H
with options: Movable on
twin castors with central
lock (1180N), paper roll
holder (1270K);
upholstery colour
64-711 ash grey



Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- In head part: Cut-out for echocardiography and upholstery part to close (in patient supine position left)
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless upwards (+50°) adjustable by gas spring
- Length head part: 800 mm
- In head part: Cut-out for echocardiography and upholstery part to close, in patient supine position left, measures: L. 300 mm x W. 250 mm
- Fixed middle part, length: 300 mm
- Foot part stepless upwards (+30°) adjustable by gas spring
- Length foot part: 850 mm
- Upholstery thickness: 60 mm



Technical data:


- Manual foot operated height adjustment: 560 mm to 980 mm by hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: 700 mm, **optional** 800 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

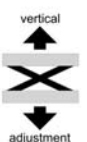

Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors \varnothing 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, \varnothing 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, \varnothing 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Paper roll holder



Art. no. 1270K / 1280K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



Art. no. 1270F / 1280F
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

Add-on elements



Art. no. 1410
1 piece of side rail, length: 350 mm



Art. no. 1850
Storage basket

Side guards



Art. no. 1080-50
1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 1060 mm, length on cut-out side: 820 mm



Art. no. 1040
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm



Art. no. 1005
1 side guard to fit on side rail, with attachment clamp, chromium-plated, length: 800 mm

Wall protection bumpers



Art. no. 1170
1 wall protection bumper \varnothing 150 mm

Upholstery options

Art. no. 2552-800
Width of upholstery: 800 mm

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

Tables for echocardiography


3-sectioned table for echocardiography with electric motorized height adjustment

... for examination, treatment ...

Model 2690-00XLE white powder coated (RAL 9010)

Model 2690-03XLE white aluminium powder coated (RAL 9006)

Model 2690-04XLE grey aluminium powder coated (RAL 9007)

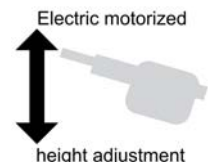

German
manufacture

All pictures:
Model 2690-04XLE
with option:
Movable on twin castors with
central lock 1187; with
special feature: foot bar;
upholstery colour
64-175 silver grey



Basic equipment:

- Table movable on twin castors, \varnothing 75 mm, each with lock
- In middle part: Cut-out for echocardiography and upholstery part to close (in patient supine position left)
- Middle part stepless positive and negative adjustable by electromotive drive, control by hand switch
- Possible positions including: Roof position, Trendelenburg, Reverse Trendelenburg
- Plug connector for potential equalization in acc. to DIN 42801
- Electromotive drive control of height adjustment by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless positive (+45°) and negative (-30°) adjustable by gas spring
- Length head part: 400 mm, with nose opening and upholstery part to close
- Middle part (length: 580 mm) stepless positive (+15°) and negative (-15°) adjustable by electromotive drive
- Middle part with cut-out for echocardiography and upholstery part to close, measures: L. 300 mm x W. 250 mm
- Foot part stepless positive (+62°) and negative (-12°) adjustable by gas spring
- Length foot part: 970 mm
- Upholstery thickness: 60 mm



Technical data:


- Electronical height adjustment: 630 mm to 1070 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 220 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandora Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

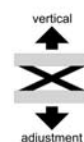

Available as an option (see pricelist for additional options)

Mobility



Art. no. 1182
Without basic mobility



Art. no. 1187
Movable on twin castors with central lock,
Ø 100 mm, lower underframe,
table height 540 mm – 980 mm



Additional motor control options



Art. no. 1920E
Hand switch (low voltage),
with magnet on the back to fasten




Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at
underframe

Nose opening

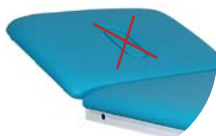

Art. no. 1615-EF
Without nose opening

Wall protection bumpers



Art. no. 1175
4 wall protection bumpers Ø 150 mm,
including mounting, for width 700 mm

Paper roll holder



Art. no. 1270F/K
Paper roll holder, mounted under head
part (not protruding),
for width 700 mm



Art. no. 1270F
Paper roll holder, mounted at foot part
(protruding: 70 mm),
for width 700 mm

Add-on elements



Art. no. 1410
1 piece of side rail, length: 350 mm



Art. no. 1005
1 side guard to fit on side rail, with
attachment clamp, chromium-plated,
length: 800 mm

Upholstery options

Art. no. 1772
Upholstery medium-soft (cold foam
cover), for width 700 mm,
thickness 70 mm

Tables for echocardiography


3-sectioned table for echocardiography with electric motorized height adjustment

... for examination, treatment ...

Model 2695-00XLE white powder coated (RAL 9010)

Model 2695-03XLE white aluminium powder coated (RAL 9006)

Model 2695-04XLE grey aluminium powder coated (RAL 9007)

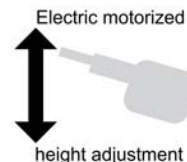

All pictures:
Model 2695-00XLE
with option: Paper roll
holder (1270F/K);
upholstery colour
64-717 deep sea



Basic equipment:

- Table movable on twin castors, \varnothing 75 mm, each with lock
- In middle part: 2 cut-outs for echocardiography and upholstery parts to close
- Middle part stepless positive and negative adjustable by electromotive drive, control by hand switch
- Possible positions including: Roof position, Trendelenburg, Reverse Trendelenburg
- Plug connector for potential equalization in acc. to DIN 42801
- Electromotive drive control of height adjustment by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless positive (+45°) and negative (-30°) adjustable by gas spring
- Length head part: 400 mm, with nose opening and upholstery part to close
- Middle part (length: 580 mm) stepless positive (+20°) and negative (-10°) adjustable by electromotive drive
- Middle part with cut-out for echocardiography and upholstery part to close, measures: L. 300 mm x W. 250 mm
- Foot part stepless positive (+60°) and negative (-12°) adjustable by gas spring
- Length foot part: 970 mm
- Upholstery thickness: 60 mm

Technical data:

- Electronical height adjustment: 630 mm to 1070 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 220 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 1182
Without basic mobility



Art. no. 1187
Movable on twin castors with central lock,
Ø 100 mm, lower underframe,
table height 540 mm – 980 mm



Additional motor control options



Art. no. 1920E
Hand switch (low voltage),
with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at
underframe

Nose opening



Art. no. 1615-EF
Without nose opening

Paper roll holder



Art. no. 1270F/K
Paper roll holder, mounted under head
part (not protruding),
for width 700 mm



Art. no. 1270F
Paper roll holder, mounted at foot part
(protruding: 70 mm),
for width 700 mm

Wall protection bumpers



Art. no. 1175
4 wall protection bumpers Ø 150 mm,
including mounting, for width 700 mm

Add-on elements



Art. no. 1410
1 piece of side rail, length: 350 mm



Art. no. 1005
1 side guard to fit on side rail, with
attachment clamp, chromium-plated,
length: 800 mm

Upholstery options

Art. no. 1772
Upholstery medium-soft
(cold foam cover),
for width 700 mm, thickness 70 mm



Tilt table for echocardiography

Tilt table with electric motorized height adjustment and tilt function
... for examination, treatment ...


Model 2095-00E white powder coated (RAL 9010)

Model 2095-03E white aluminium powder coated (RAL 9006)

Model 2095-04E grey aluminium powder coated (RAL 9007)



All pictures:
Model 2095-00E
with options:
1 pair of belt guides
(2090-036), 1 safety belt with
automatic lock (2090-0351),
upholstery for strap (2090-034);
upholstery colour: 64-702 white

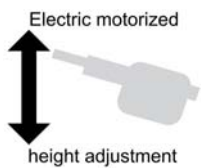

Basic equipment:

- Lying surface 2-sectioned, cut-out for echocardiography with upholstery part to close in head part
- Operation of electric motorized height adjustment by foot switch, operation of electric motorized tilt function by hand switch (independent from each other)
- Integrated safety function for automatic blocking of electromotive drive
- Height adjustment by electromotive drive with freewheeling downwards, automatic declutch in case of resistance
- Plug connector for potential equalization in acc. to DIN 42801
- Movable on 4 twin castors Ø 75 mm, each with lock
- Foot plate to stand on (upholstered), removable
- Side rails along lying surface (left and right side) for special tilt table equipment
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless adjustable from 0° to +45° by gas spring, length head part: 800 mm
- In head part: Cut-out for echocardiography, measures: Length 300 mm x Width 250 mm
- Cut-out for echocardiography and upholstery part to close (in patient supine position left)
- Fixed foot part, length foot part: 1130 mm
- Thickness of upholstery: 80 mm

Technical data:


- Electric motorized height adjustment: 620 mm to 930 mm in 34 sec. by AC-Motor, pressure capacity of 10.000 N
- Electric motorized tilt function from 0° to +85° in 22 sec. by AC-Motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm | Width overall: 850 mm
- Length of lying surface: 1930 mm | Length overall: 1960 mm
- Load capacity: 180 kg; according to IEC 60601 tested with fourfold static load

kg

Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

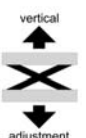

Attention: necessary room height: 2.45 m

Available as an option (see pricelist for additional options)

Additional motor control options



Art. no. 1920E
Hand switch (low voltage),
with magnet on the back to
fasten



Art. no. 1910E
1 foot switch (low voltage)

Grip handles



Art. no. 2090-028
1 pair grip handles, small




Art. no. 2090-038
1 pair grip handles, big

Belt elements



Art. no. 2090-034
Upholstery for strap

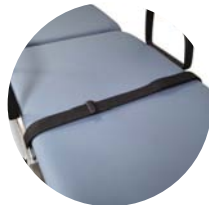

Art. no. 2090-035
1 safety belt with velcro fastener



Art. no. 2090-0351
1 safety belt with automatic
lock



Art. no. 2090-036
1 pair belt guides

Arm rest



Art. no. 2090-0143
Attachment clamp to fix
arm rest at side rail

Art. no. 2090-059
Arm rest

Paper roll holder




Art. no. 1272
Paper roll holder, mounted under head
part (not protruding)

Special equipment



Art. no. 2090-090
Angle gauge, on both sides
of the table

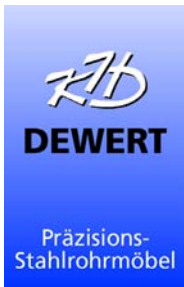

K.H. Dewert GmbH
Vollmestraße 7
33649 Bielefeld
Germany

Post office 14 40 28
33634 Bielefeld
Germany

Phone + 49 (0) 5 21 / 4 00 27-0
Fax + 49 (0) 5 21 / 4 00 27-27

Internet: www.khdewert.de
E-Mail: info@khdewert.de
YouTube: K.H. Dewert

TABLES FOR ECHOCARDIOGRAPHY FOR EXAMINATION AND TREATMENT

TABLES FOR SURGERY AND HOSPITAL

