


VERSATILE TABLES
FOR EXAMINATION, TREATMENT AND THERAPY



SERIES 2000XL AND 4000E

2023.1

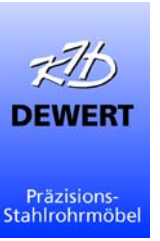

VERSATILE TABLES ...



... WITH
INCLINATION ADJUSTMENT




... WITH
VARIABLE LENGTH



SHORT TABLES



TABLES FOR SURGERY AND HOSPITAL

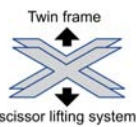

Series 2000XL and 4000E - Versatile tables



The models of series 2000XL and 4000E ensure optimal support in examination and treatment without limitation.

A proven and reliable height adjustment is either effected by a robust hydraulic pump with an optimal ratio or a high-performance freewheeling electric motor with high adjustment speed which **corresponds to the newest advice of BfArM.**

In addition, Series 4000E with electromotive drive has a high adjustment range. Due to the extremely stable twin frame scissor lifting system a high load capacity of up to 280 kg (depending on model) is possible.

Equipped with our proven upholstery made of laminated wood with homogenous polyurethane foam, the models of versatile tables enable a pleasant resting comfort.

Furthermore, every table is equipped with the Dewert-hygiene-package, which guarantees optimal cleaning and disinfection of the upholstery as well as of the frame construction.

Therefore, the high quality-models of series 2000XL and 4000E have best prerequisites for long-lasting daily practice use - without compromise.




Content



Good reasons for a versatile table	2
Versatile tables for examination, treatment and therapy	
2-sectioned versatile tables:	
• Model 2000XLE / 2001XLE	6
• Model 2000XL/H / 2001XL/H	8
• Model 4000XLE / 4001XLE	10
• Model 2030XLE / 2035XLE, long head part	12
• Model 2030XL/H / 2035XL/H, long head part	14
2-sectioned versatile tables with inclination adjustment:	
• Model 2052XLE / 2053XLE	16
• Model 2052XL/H / 2053XL/H	18
• Model 2054XL/H, long head part	20
Versatile tables for examination, treatment and therapy with variable length	
2-sectioned versatile tables with variable length:	
• Model 2006XLE / 2007XLE	22
• Model 2006XL/H / 2007XL/H	24
3-sectioned versatile tables with variable length:	
• Model 2008XLE / 2009XLE	26
• Model 2008XL/H / 2009XL/H	28
• Model 2010XLE / 2011 XLE	30
• Model 2010XL/H / 2011XL/H	32
Versatile tables for examination, treatment and therapy Short tables	
• Model 2100XLE / 2105XLE / 2110XLE	34
• Model 2100XL/H / 2105XL/H / 2110XL/H	36
• Model 2114XLE / 2115XLE	38
• Model 2114XL/H / 2115XL/H	40


Series 2000XL and 4000E - Versatile tables

Good reasons for a versatile table

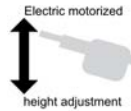
Electric motorized height adjustment



- Electronical height adjustment by AC-motor with a pressure capacity of 10000 N
- Constant adjustment speed even in case of high load
- EPR-safety power cable with strain relief, length: 4 m
- **Integrated safety devices (according to the newest advice of BfArM)**
 - **Automatic blocking** of electromotive drive 
 - **Active freewheeling downwards**, automatic declutch in case of resistance 
- Series 4000E: high adjustment range
- Load capacity: 250 kg; Model 2000XLE: 280 kg

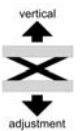

Hydraulic height adjustment

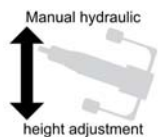

- Manual foot operated height adjustment by a hydraulic pump with a pressure capacity of 10000 N
- Pedals for foot-controlled height positioning on **both sides** of the table
- Load capacity: 250 kg



Dewert-hygiene-equipment


Enables optimal cleaning and disinfection



- Hinge lid made of identical upholstery covering
- Artificial leather underneath the upholstery
- Vent holes underneath the upholstery, closed by special caps with valve function

- Easy hygiene because of an open construction: good accessibility of components, no plastic coverings which are breakable in the ageing process

- Option: Innovative artificial leather Skai® Toronto EN with staynu

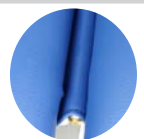

Series 2000XL and 4000E - Versatile tables

Good reasons for a versatile table




Durable ✓



Robust ✓



Easy-care ✓

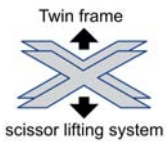


Twin frame scissor lifting system ✓



High load capacity ✓



Proven Dewert-quality

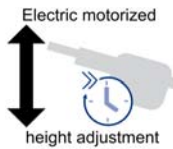

High adjustment speed ✓



Universal use ✓

Technology:



- Extremely stable twin frame scissor lifting system, with high torsional stiffness due to thick-walled square tubes with big dimensions, also for under and upper frame and scissor

- High stability due to massive 20 mm axes




- Widely-spaced feet and a high deadweight guarantee stability

- Robust frame construction which is proven in demanding everyday life of hospitals for years



- Abrasion-resistant guide system at under and upper frame due to the use of roller bearings

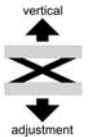

Upholstery:

- Thickness of upholstery: 60 mm
- Multiple glued, high-strength laminated wood panel (thickness: 18 mm) no breaking chipboards in case of pressure on the edges



- Two-layer homogenous polyurethane foam, density 46 and 75, Ökotex-certified

- 2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Cleaning kit for Skai products

(Art. no. 050-1000)

Consists of:

- Spray bottle with 200 ml special Skai cleansing agent
- 1 microfibre care cloth 40 x 40 cm

Series 2000XL and 4000E - Versatile tables

Good reasons for a versatile table

... with 2- or 3-sectioned lying surface:



- Adjustable head part either positive and negative or only positive adjustable with tubular steel frame for high torsion rigidity
- Pressure stable gas springs ensure working at head part and foot part **without vibrations**
- Operating lever to manipulate the gas spring easily accessible
- 3-sectioned lying surface in addition:
Variable length of table due to foot part which can be hinged vertical



... with 2-sectioned lying surface and inclination adjustment:



- Lying surface inclinable by gas spring or electromotive drive (depending on model)
- **Models 2052XLE and 2053XLE:**
 - Lying surface inclinable to 22° (Trendelenburg) by electromotive drive in 18 seconds
 - Operation by hand switch



- **Models 2052XL/H, 2053XL/H and 2054XL/H:**
 - Lying surface inclinable to 12° (Trendelenburg) by gas spring
 - A push handle at foot part facilitates the operation of manual inclination



... with 1-sectioned lying surface:



- **Short table**, with 3-side upholstered border (depending on model)

Series 2000XL


2-sectioned versatile treatment tables with electric motorized height adjustment
... for examination, treatment, therapy ...

Model 2000-00XLE and 2001-00XLE white powder coated (RAL 9010)

Model 2000-03XLE and 2001-03XLE white aluminium powder coated (RAL 9006)

Model 2000-04XLE and 2001-04XLE silver aluminium powder coated (RAL 9007)

► Model 2805-..XLE corresponds with model 2000-..XLE, but with nose opening

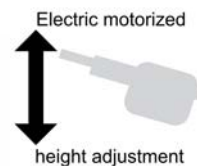

Model 2000-00XLE
with options: Movable on
twin castors with central lock
(1180N), paper roll holder
(1270K);
upholstery colour: 64-702 white



Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive (art. no. 1999, mandatory equipment)
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless up- (+40°) and downwards (-35°) adjustable by gas spring
- Length head part: 550 mm
- Fixed foot part
- Length foot part: 1400 mm
- Upholstery thickness: 60 mm



Model 2000-04XLE with options:
Wheel-lifting/fixing mechanism (1045N), foot switch rail, surround
(1960E), nose opening and upholstery part to close (1610, 1620)

Technical data:

- Electronical height adjustment: 480 mm to 920 mm in 22 sec.
by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: Model 2000-..XLE: 700 mm | Model 2001-..XLE: 800 mm
- Length of lying surface: 1950 mm
- Load capacity: 280 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandora Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Additional motor control options



Art. no. 1920E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at underframe



Art. no. 1950E
Foot switch rail, lateral, on both sides, security adjustment



Art. no. 1960E
Foot switch rail, surround, security adjustment

Paper roll holder



Art. no. 1270K / 1280K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



Art. no. 1270F / 1280F
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close

Side guards



Art. no. 1080
1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 1060 mm



Art. no. 1040
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm

Wall protection bumpers



Art. no. 1175 / 1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Upholstery options

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

Series 2000XL


2-sectioned versatile treatment tables with hydraulic height adjustment
... for examination, treatment, therapy ...

Model 2000-00XL/H and 2001-00XL/H white powder coated (RAL 9010)

Model 2000-03XL/H and 2001-03XL/H white aluminium powder coated (RAL 9006)

Model 2000-04XL/H and 2001-04XL/H silver aluminium powder coated (RAL 9007)

▶ **Model 2805-..XL/H corresponds with model 2000-..XL/H, but with nose opening**

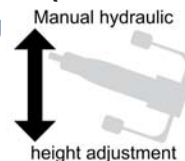

Both pictures:
Model 2000-00XL/H
with options:
Movable on twin castors
with central lock (1180N),
push handle (picture
shows chromium-plated
push handle, basic push
handle is RAL 9006 pow-
der coated);
upholstery colour:
64-722 orange

Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless up- (+40°) and downwards (-35°) adjustable by gas spring
- Length head part: 550 mm
- Fixed foot part
- Length foot part: 1400 mm
- Upholstery thickness: 60 mm



Technical data:

- Manual foot operated height adjustment: 500 mm to 920 mm by a hydraulic pump with a pressure capacity of 10000 N
- Width of lying surface: Model 2000-..XL/H: 700 mm | 2001-..XL/H: 800 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandora Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Paper roll holder



Art. no. 1270K / 1280K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



Art. no. 1270F / 1280F
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

Wall protection bumpers



Art. no. 1175 / 1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Side guards



Art. no. 1080
1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 1060 mm



Art. no. 1040
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close

Upholstery options

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

Series 4000E

2-sectioned versatile treatment tables with electric motorized height adjustment
... for examination, treatment, therapy ...

Model 4000-00E and 4001-00E white powder coated (RAL 9010)

Model 4000-03E and 4001-03E white aluminium powder coated (RAL 9006)

Model 4000-04E and 4001-04E silver aluminium powder coated (RAL 9007)

► Model 4805-..E corresponds with model 4000-..E, but with nose opening

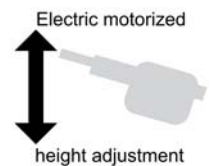


Table in highest position

All pictures:
Model 4000-00XLE
with options:
Wheel-lifting/fixing mechanism
(1045N), without foot switch
(1910E-EF), hand switch (1920E);
upholstery colour: 64-173 smaragd

Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- High adjustment range
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive (art. no. 1999, mandatory equipment)
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless up- (+40°) and downwards (-35°) adjustable by gas spring
- Length head part: 500 mm
- Fixed foot part
- Length foot part: 1500 mm
- Upholstery thickness: 60 mm



Table in lowest position

Technical data:


- Electronical height adjustment: 420 mm to 970 mm in 18 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: Model 4000E: 700 mm | Model 4001E: 800 mm
- Length of lying surface: 2000 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

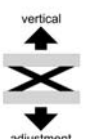

Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Additional motor control options



Art. no. 1920E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at underframe



Art. no. 1950E
Foot switch rail, lateral, on both sides, security adjustment



Art. no. 1960E
Foot switch rail, surround, security adjustment

Paper roll holder



Art. no. 1270K / 1280K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



Art. no. 1270F / 1280F
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close

Side guards



Art. no. 1080
1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 1060 mm

Wall protection bumpers



Art. no. 1175 / 1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Upholstery options

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm


Series 2000XL

2-sectioned versatile treatment tables with electric motorized height adjustment
... for examination, treatment, therapy ...

Model 2030-00XLE and 2035-00XLE white powder coated (RAL 9010)

Model 2030-03XLE and 2035-03XLE white aluminium powder coated (RAL 9006)

Model 2030-04XLE and 2035-04XLE silver aluminium powder coated (RAL 9007)

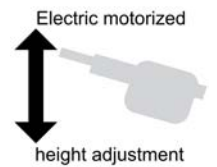

Both pictures:
Model 2030-04XLE
with options: Movable on
twin castors with central lock
(1180N), paper roll holder
(1270F);
upholstery colour:
64-717 deep sea



Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Long head part with half-round design
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive (art. no. 1999, mandatory equipment)
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless adjustable from -20° to +70° by gas spring (sitting position)
- Length head part: 750 mm
- Fixed foot part
- Length foot part: 1.200 mm
- Upholstery thickness: 60 mm



Technical data:


- Electronical height adjustment: 480 mm to 920 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: Model 2030-...XLE: 700 mm | Model 2035-...XLE: 800 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

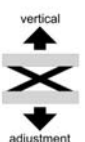

Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering

Additional motor control options



Art. no. 1920E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at underframe



Art. no. 1950E
Foot switch rail, lateral, on both sides, security adjustment

Paper roll holder



Art. no. 1270K / 1280K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



Art. no. 1270F / 1280F
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close

Side guards



Art. no. 1080-60
1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 820 mm

Wall protection bumpers



Art. no. 1175 / 1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Upholstery options

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

Series 2000XL


2-sectioned versatile treatment tables with hydraulic height adjustment

... for examination, treatment, therapy ...

Model 2030-00XL/H and 2035-00XL/H white powder coated (RAL 9010)


Model 2030-03XL/H and 2035-03XL/H white aluminium powder coated (RAL 9006)

Model 2030-04XL/H and 2035-04XL/H silver aluminium powder coated (RAL 9007)

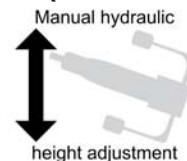

All pictures:
Model 2035-00XL/H
with options: Movable on
twin castors with central
lock (1180N), paper roll
holder (1280K);
upholstery colour:
64-717 deep sea

Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Long head part with half-round design
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment 


Upholstery:

Lying surface 2-sectioned:

- Head part stepless adjustable from -20° to $+70^\circ$ by gas spring (sitting position)
- Length head part: 750 mm
- Fixed foot part
- Length foot part: 1.200 mm
- Upholstery thickness: 60 mm



Technical data:

- Manual foot operated height adjustment: 540 mm to 960 mm by a hydraulic pump with a pressure capacity of 10000 N
- Width of lying surface: Model 2030-..XL/H: 700 mm | Model 2035-..XL/H: 800 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load 

Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Paper roll holder



Art. no. 1270K / 1280K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



Art. no. 1270F / 1280F
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

Wall protection bumpers



Art. no. 1175 / 1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Side guards



Art. no. 1080-60
1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 820 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close

Upholstery options

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

Series 2000XL


2-sectioned versatile treatment tables with electric motorized height adjustment
... for examination, treatment, therapy ...

Model 2052-00XLE and 2053-00XLE white powder coated (RAL 9010)

Model 2052-03XLE and 2053-03XLE white aluminium powder coated (RAL 9006)

Model 2052-04XLE and 2053-04XLE silver aluminium powder coated (RAL 9007)

► Model 2910-..XLE corresponds with model 2052-..XLE, but with nose opening

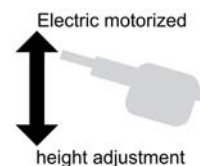

3 upper pictures:
Model 2052-00XLE
with options:
Movable on comfort twin
castors with central lock
(1185); upholstery colour:
64-716 black

Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Lying surface inclinable to 22° (Trendelenburg) by electromotive drive in 18 sec.
- Operation of electric motor to adjust lying surface by 1 hand switch
- Operation of electric motor for height adjustment by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drives (art. no. 1999, mandatory equipment)
- Electric motor for height adjustment: Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless up- (+45°) and downwards (-20°) adjustable by gas spring
- Length head part: 500 mm, optional with nose opening and upholstery part
- Fixed foot part
- Length foot part: 1.500 mm
- Upholstery thickness: 60 mm



Model 2052-00XLE with options: Movable on twin castors with central lock (1180N), 2 side guards lateral lowerable (1080); upholstery colour: 64-174 sky

Technical data:


- Electronical height adjustment: 540 mm to 980 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Lying surface inclinable to 22° (Trendelenburg) by electromotive drive in 18 sec.
- Width of lying surface: Model 2052-..XLE: 700 mm | Model 2053-..XLE: 800 mm
- Length of lying surface: 2000 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandora Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

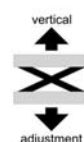

Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Additional motor control options



Art. no. 1920E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at underframe



Art. no. 1950E
Foot switch rail, lateral, on both sides, security adjustment

Paper roll holder



Art. no. 1270K / 1280K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



Art. no. 1270F / 1280F
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close

Side guards



Art. no. 1080
1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 1060 mm

Upholstery options

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

Series 2000XL


2-sectioned versatile treatment tables with hydraulic height adjustment
... for examination, treatment, therapy ...

Model 2052-00XL/H and 2053-00XL/H white powder coated (RAL 9010)

Model 2052-03XL/H and 2053-03XL/H white aluminium powder coated (RAL 9006)

Model 2052-04XL/H and 2053-04XL/H silver aluminium powder coated (RAL 9007)

▶ **Model 2810-..XL/H corresponds with model 2052-..XL/H, but with nose opening**

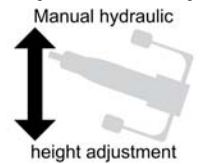

All pictures:
Model 2052-04XL/H
with option: Movable on com-
fort twin castors with central
lock (1185) paper roll holder
(1270K), 2 side guards lateral
lowerable (1080) (picture
shows chromium-plated push
handle, basic push handle is
RAL 9006 powder coated);
upholstery colour:
64-717 deep sea

Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Lying surface inclinable to 12° (Trendelenburg) by gas spring
- Push handle at foot part facilitates the operation of manual inclination, RAL 9006 powder coated
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless up- (+45°) and downwards (-35°) adjustable by gas spring
- Length head part: 500 mm
- Fixed foot part
- Length foot part: 1500 mm
- Upholstery thickness: 60 mm



Technical data:


- Manual foot operated height adjustment: 550 mm to 970 mm by a hydraulic pump with a pressure capacity of 10000 N
- Lying surface inclinable to 12° (Trendelenburg) by gas spring
- Width of lying surface: Model 2052-..XL/H: 700 mm | Model 2053-..XL/H: 800 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandora Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

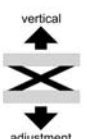

Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Paper roll holder



Art. no. 1270K / 1280K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm

Wall protection bumpers



Art. no. 1175 / 1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close

Side guards



Art. no. 1080
1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 1060 mm

Upholstery options

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for width 700 mm/ 800 mm, thickness 80 mm

Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm


Series 2000XL

2-sectioned versatile treatment tables with hydraulic height adjustment
... for examination, treatment, therapy ...

Model 2054-00XL/H white powder coated (RAL 9010)

Model 2054-03XL/H white aluminium powder coated (RAL 9006)

Model 2054-04XL/H silver aluminium powder coated (RAL 9007)

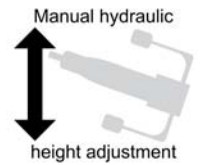

All pictures:
Model 2054-00XL/H
with option: Movable on
twin castors with central
lock (1180N) (picture
shows chromium-plated
push handle, basic
push handle is RAL
9006 powder coated);
upholstery colour:
64-703 pearl

Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Lying surface inclinable to 12° (Trendelenburg) by gas spring
- Push handle at foot part facilitates the operation of manual inclination, RAL 9006 powder coated
- Long head part with half-round design
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless adjustable from 0° to +75° by gas spring (sitting position)
- Length head part: 750 mm
- Fixed foot part
- Length foot part: 1200 mm
- Upholstery thickness: 60 mm



Technical data:


- Manual foot operated height adjustment: 550 mm to 970 mm by a hydraulic pump with a pressure capacity of 10000 N
- Lying surface inclinable to 12° (Trendelenburg) by gas spring
- Width of lying surface: 700 mm, optional 800 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

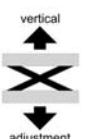

Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Paper roll holder



Art. no. 1270F/K / 1280F/K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm

Wall protection bumpers



Art. no. 1175 / 1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close

Side guards



Art. no. 1080
1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 1060 mm

Upholstery options

Art. no. 2054-800
Width of upholstery 800 mm

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

Series 2000XL


2-sectioned versatile treatment tables with electric motorized height adjustment
... for examination, treatment, therapy ...

Model 2006-00XLE and 2007-00XLE white powder coated (RAL 9010)

Model 2006-03XLE and 2007-03XLE white aluminium powder coated (RAL 9006)

Model 2006-04XLE and 2007-04XLE silver aluminium powder coated (RAL 9007)

► Model 2806-..XLE corresponds with model 2006-..XLE, but with nose opening

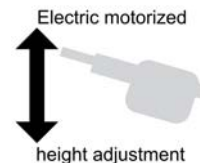

All pictures:
Model 2006-00XLE
with option: Movable on
comfort twin castors,
antistatic (1186-10);
upholstery colour:
64-174 sky

Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Head part can be hinged vertical
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive (art. no. 1999, mandatory equipment)
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless up- (+30°) and downwards (-90°) adjustable by gas spring
- Length head part: 450 mm
- Fixed foot part
- Length foot part: 1500 mm
- Upholstery thickness: 60 mm



Technical data:


- Electronical height adjustment: 520 mm to 960 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: Model 2006-..XLE: 700 mm | Model 2007-..XLE: 800 mm
- Length of lying surface: min. 1550 mm - max. 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

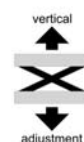

Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Additional motor control options



Art. no. 1920E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at underframe



Art. no. 1950E
Foot switch rail, lateral, on both sides, security adjustment

Paper roll holder



Art. no. 1270K / 1280K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



Art. no. 1270F / 1280F
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

Side guards



Art. no. 1080
1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 1060 mm



Art. no. 1040
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm

Wall protection bumpers



Art. no. 1175 / 1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close

Upholstery options

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

Series 2000XL


2-sectioned versatile treatment tables with hydraulic height adjustment
... for examination, treatment, therapy ...

Model 2006-00XL/H and 2007-00XL/H white powder coated (RAL 9010)

Model 2006-03XL/H and 2007-03XL/H white aluminium powder coated (RAL 9006)

Model 2006-04XL/H and 2007-04XL/H silver aluminium powder coated (RAL 9007)

▶ **Model 2806-..XL/H corresponds with model 2006-..XL/H, but with nose opening**

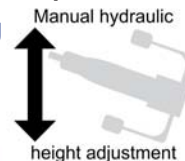

All pictures:
Model 2007-00XL/H
with option: Paper roll
holder (1280K);
upholstery colour:
64-721 atoll



Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 2-sectioned:

- Head part stepless up- (+30°) and downwards (-90°) adjustable by gas spring
- Length head part: 450 mm
- Fixed foot part
- Length foot part: 1500 mm
- Upholstery thickness: 60 mm



Technical data:

- Manual foot operated height adjustment: 540 mm to 960 mm by a hydraulic pump with a pressure capacity of 10000 N
- Width of lying surface: Model 2006-..XL/H: 700 mm | Model 2007-..XL/H: 800 mm
- Length of lying surface: min. 1550 mm - max. 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandora Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Paper roll holder



Art. no. 1270K / 1280K
Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



Art. no. 1270F / 1280F
Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

Wall protection bumpers



Art. no. 1175 / 1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Side guards



Art. no. 1080
1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 1060 mm



Art. no. 1040
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close

Upholstery options

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

Series 2000XL


3-sectioned versatile treatment tables with electric motorized height adjustment ... for examination, treatment, therapy ...

Model 2008-00XLE and 2009-00XLE white powder coated (RAL 9010)

Model 2008-03XLE and 2009-03XLE white aluminium powder coated (RAL 9006)

Model 2008-04XLE and 2009-04XLE silver aluminium powder coated (RAL 9007)

▶ Model 2808-..XLE corresponds with model 2008-..XLE, but with nose opening

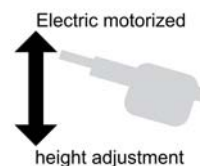

Both pictures:
Model 2008-00XLE
with options:
Movable on twin castors with central lock (1180N), nose opening and upholstery part to close (1610, 1620);
upholstery colour:
64-701 stone blue

Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Foot part can be hinged vertical
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive (art. no. 1999, mandatory equipment)
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+50°) and downwards (-45°) adjustable by gas spring
- Length head part: 450 mm
- Fixed middle part
- Length middle part: 1290 mm
- Foot part stepless up- (+30°) and downwards (-90°) adjustable by gas spring
- Length foot part: 450 mm
- Upholstery thickness: 60 mm



Technical data:


- Electronical height adjustment: 520 mm to 960 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: Model 2008-..XLE: 700 mm | Model 2009-..XLE: 800 mm
- Length of lying surface: min. 1790 mm - max. 2190 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

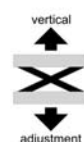

Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Additional motor control options



Art. no. 1920E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at underframe



Art. no. 1950E
Foot switch rail, lateral, on both sides, security adjustment

Paper roll holder



Art. no. 1270K/F / 1280K/F
Paper roll holder, mounted at foot part (not protruding), for width 700 mm / 800 mm

Wall protection bumpers



Art. no. 1175 / 1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close

Upholstery options

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

Series 2000XL


3-sectioned versatile treatment tables with hydraulic height adjustment
... for examination, treatment, therapy ...

Model 2008-00XL/H and 2009-00XL/H white powder coated (RAL 9010)

Model 2008-03XL/H and 2009-03XL/H white aluminium powder coated (RAL 9006)

Model 2008-04XL/H and 2009-04XL/H silver aluminium powder coated (RAL 9007)

► **Model 2808-..XL/H corresponds with model 2008-..XL/H, but with nose opening**

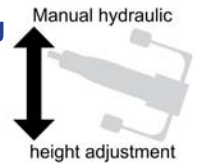

All pictures:
Model 2008-00XL/H
with option:
Nose opening and
upholstery part to close
(1610 and 1620);
upholstery colour: 64-717 deep sea



Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Foot part can be hinged vertical
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+50°) and downwards (-45°) adjustable by gas spring
- Length head part: 450 mm
- Fixed middle part
- Length middle part: 1290 mm
- Foot part stepless up- (+30°) and downwards (-90°) adjustable by gas spring
- Length foot part: 450 mm
- Upholstery thickness: 60 mm



Technical data:


- Manual foot operated height adjustment: 540 mm to 960 mm by a hydraulic pump with a pressure capacity of 10000 N
- Width of lying surface: Model 2008-..XL/H: 700 mm | Model 2009-..XL/H: 800 mm
- Length of lying surface: min. 1790 mm - max. 2190 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandora Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

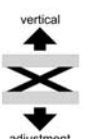

Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Paper roll holder



Art. no. 1270K/F / 1280K/F
Paper roll holder, mounted at foot part (not protruding), for width 700 mm / 800 mm

Wall protection bumpers



Art. no. 1175 / 1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close

Upholstery options

Art. no. 1770 / 1780

Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782

Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

Series 2000XL


3-sectioned versatile treatment tables with with electric motorized height adjustment ... for examination, treatment, therapy ...

Model 2010-00XLE and 2011-00XLE white powder coated (RAL 9010)

Model 2010-03XLE and 2011-03XLE white aluminium powder coated (RAL 9006)

Model 2010-04XLE and 2011-04XLE silver aluminium powder coated (RAL 9007)

► **Model 2807-..XLE corresponds with model 2010-..XLE, but with nose opening**

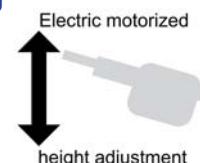

Model 2010-00XLE
with options:
Movable on comfort twin castors
with central lock (1185);
upholstery colour:
64-717 ddep sea

Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Foot part can be hinged vertical
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive (art. no. 1999, mandatory equipment)
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless adjustable from 0° to +60° by gas spring
- Length head part: 450 mm
- Fixed middle part
- Length middle part: 1.050 mm
- Foot part stepless up- (+30°) and downwards (-90°) adjustable by gas spring
- Length foot part: 450 mm
- Upholstery thickness: 60 mm



Model 2010-00XLE with options:
Movable on comfort twin castors with central lock (1185), paper roll holder (1270K/F)

Technical data:

- Electronical height adjustment: 530 mm to 970 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: Model 2010-..XLE: 700 mm | Model 2011-..XLE: 800 mm
- GesamtLength of lying surface: min. 1550 mm - max. 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Additional motor control options



Art. no. 1920E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at underframe



Art. no. 1950E
Foot switch rail, lateral, on both sides, security adjustment

Paper roll holder



Art. no. 1270F/K / 1280F/K
Paper roll holder, mounted under head part (70 mm protruding), for width 700 mm / 800 mm



Art. no. 1270K/F / 1280K/F
Paper roll holder, mounted at foot part (not protruding), for width 700 mm / 800 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close

Wall protection bumpers



Art. no. 1175 / 1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Side guards



Art. no. 1080-60
1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 820 mm

Upholstery options

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

Series 2000XL


3-sectioned versatile treatment tables with hydraulic height adjustment
... for examination, treatment, therapy ...

Model 2010-00XL/H and 2011-00XL/H white powder coated (RAL 9010)

Model 2010-03XL/H and 2011-03XL/H white aluminium powder coated (RAL 9006)

Model 2010-04XL/H and 2011-04XL/H silver aluminium powder coated (RAL 9007)


▶ Model 2807-..XL/H corresponds with model 2010-..XL/H, but with nose opening

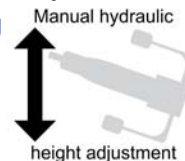

Both pictures:
Model 2010-00XL/H
with option: Movable on twin
castors with central lock (1180N),
paper roll holder (1270K/F);
upholstery colour: 64-708 cherry



Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Foot part can be hinged vertical
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment 

Upholstery:

Lying surface 3-sectioned:

- Head part stepless adjustable from 0° to +60° by gas spring
- Length head part: 450 mm
- Fixed middle part
- Length middle part: 1050 mm
- Foot part stepless up- (+30°) and downwards (-90°) adjustable by gas spring
- Length foot part: 450 mm
- Upholstery thickness: 60 mm



Technical data:

- Manual foot operated height adjustment: 550 mm to 970 mm by a hydraulic pump with a pressure capacity of 10000 N
- Width of lying surface: Model 2010-..XL/H: 700 mm | Model 2011-..XL/H: 800 mm
- Length of lying surface: min. 1550 mm - max. 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandora Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Paper roll holder



Art. no. 1270F/K / 1280F/K
Paper roll holder, mounted under head part (70 mm protruding), for width 700 mm / 800 mm



Art. no. 1270K/F / 1280K/F
Paper roll holder, mounted at foot part (not protruding), for width 700 mm / 800 mm

Nose opening



Art. no. 1610
Nose opening, in head part

Art. no. 1620
Upholstery part to close

Wall protection bumpers



Art. no. 1175 / 1176
4 wall protection bumpers Ø 150 mm, including mounting, for width 700 mm / 800 mm

Side guards



Art. no. 1080-60
1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 820 mm

Upholstery options

Art. no. 1770 / 1780
Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

Series 2000XL


Short tables with electric motorized height adjustment

... for examination, treatment, therapy ...

Model 2100-00XLE, 2105-00XLE, 2110-00XLE white powder coated (RAL 9010)

Model 2100-03XLE, 2105-03XLE, 2110-03XLE white aluminium powder coated (RAL 9006)

Model 2100-04XLE, 2105-04XLE, 2110-04XLE silver aluminium powder coated (RAL 9007)

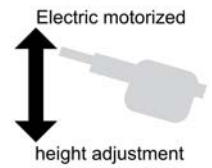

All pictures:
Model 2100-00XLE
with option: Movable on comfort twin castors, each with lock (1120N-10), without foot switch (1910E-EF), hand switch (1920E); upholstery colour: 64-717 deep sea



Table in lowest position

Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive (art. no. 1999, mandatory equipment)
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 1-sectioned:

- Length: 1400 mm (other measures on request)
- Upholstery thickness: 60 mm



Table in highest position

Technical data:


- Electronical height adjustment: 480 mm to 920 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: Model 2100-..XLE: 650 mm | Model 2105-..XLE: 800 mm | Model 2110-..XLE: 1.000 mm
- Length of lying surface: 1400 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

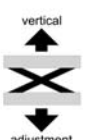

Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Additional motor control options



Art. no. 1920E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at underframe

Series 2000XL


Short tables with hydraulic height adjustment

... for examination, treatment, therapy ...

Model 2100-00XL/H, 2105-00XL/H, 2110-00XL/H white powder coated (RAL 9010)

Model 2100-03XL/H, 2105-03XL/H, 2110-03XL/H white aluminium powder coated (RAL 9006)

Model 2100-04XL/H, 2105-04XL/H, 2110-04XL/H silver aluminium powder coated (RAL 9007)

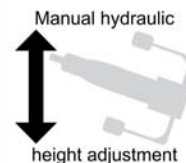

Both pictures:
Model 2100-00XL/H
with option: Movable on
twin castors with central
lock (1180N);
upholstery colour: 64-717 deep sea



Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 1-sectioned:

- Length: 1400 mm (other measures on request)
- Upholstery thickness: 60 mm



Table in highest position



Technical data:

- Manual foot operated height adjustment: 50 mm to 920 mm by a hydraulic pump with a pressure capacity of 10000 N
- Width of lying surface: Model 2100-..XL/H: 650 mm | Model 2105-..XL/H: 800 mm | Model 2110-..XL/H: 1.000 mm
- GesamtLength of lying surface: 1400 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandora Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin
castors \varnothing 100 mm,
each with total lock,
basic height -15 mm



Art. no. 1180N
Movable on twin castors,
 \varnothing 100 mm with central
lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin
castors, \varnothing 125 mm with
central lock,
basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering

Series 2000XL


Short tables with electric motorized height adjustment

... for examination, treatment, therapy ...

Model 2114-00XLE and 2115-00XLE white powder coated (RAL 9010)

Model 2114-03XLE and 2115-03XLE white aluminium powder coated (RAL 9006)

Model 2114-04XLE and 2115-04XLE silver aluminium powder coated (RAL 9007)

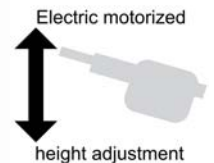

Model 2114-04XLE
with basic equipment;
upholstery colour:
64-714 olive green

Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive (art. no. 1999, mandatory equipment)
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 1-sectioned:

- Length: 1400 mm (other measures on request)
- with 3-side upholstered border, height: 150 mm
- Upholstery thickness: 60 mm



Technical data:


- Electronical height adjustment: 480 mm to 920 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: Model 2114-...XLE: 650 mm | Model 2115-...XLE: 800 mm
- GesamtLength of lying surface: 1400 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

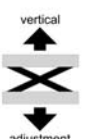

Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Additional motor control options



Art. no. 1920E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
Foot switch (low voltage)



Art. no. 1940E
1 foot switch fixation mounted at underframe

Series 2000XL


Short tables with hydraulic height adjustment

... for examination, treatment, therapy ...

Model 2114-00XL/H and 2115-00XL/H white powder coated (RAL 9010)

Model 2114-03XL/H and 2115-03XL/H white aluminium powder coated (RAL 9006)

Model 2114-04XL/H and 2115-04XL/H silver aluminium powder coated (RAL 9007)

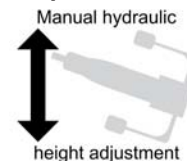

Model 2114-00XL/H
with option: Movable on
twin castors each with
lock (1120N-10);
upholstery colour: 64-717 deep sea



Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 1-sectioned:

- Length: 1400 mm (other measures on request)
- with 3-side upholstered border, height: 150 mm
- Upholstery thickness: 60 mm



Technical data:

- Manual foot operated height adjustment: 500 mm to 920 mm by a hydraulic pump with a pressure capacity of 10000 N
- Width of lying surface: Model 2114-..XL/H: 650 mm | Model 2115-..XL/H: 800 mm
- Length of lying surface: 1400 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandora Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
Movable on comfort twin
castors Ø 100 mm,
each with total lock,
basic height -15 mm




Art. no. 1180N
Movable on twin castors,
Ø 100 mm with central
lock, basic height +20 mm



Art. no. 1185
Movable on comfort twin
castors, Ø 125 mm with
central lock,
basic height + 50 mm



Art. no. 1045N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering

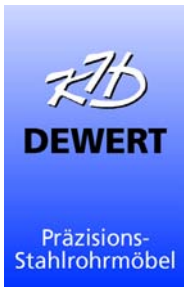

K.H. Dewert GmbH
Vollmestraße 7
33649 Bielefeld
Germany

Post office 14 40 28
33634 Bielefeld
Germany

Phone + 49 (0) 5 21 / 4 00 27-0
Fax + 49 (0) 5 21 / 4 00 27-27

Internet: www.khdewert.de
E-Mail: info@khdewert.de
YouTube: K.H. Dewert

VERSATILE TABLES FOR EXAMINATION, TREATMENT AND THERAPY

You will find the whole Dewert product range with additional articles and equipment on our website and in our further catalogues. We look forward to hear from you!



04.2023_01

TABLES FOR SURGERY AND HOSPITAL

Technical data subject to change and subject to errors

